

MICRONIC | SAMPLE STORAGE | RACKS

LASER ETCHED 1D WHITE STORAGE RACK

DO YOU REQUIRE A STORAGE RACK THAT HAS A PERMANENT AND UNIQUE BARCODE? IF SO, THEN MICRONIC OFFERS THE LASER ETCHED LOBORACK-96 AND ROBORACK-96, WHICH CAN HOLD 96 INDIVIDUAL TUBES IN VOLUMES OF 0.50ML, 0.75ML, 1.40ML AND 2.50ML IN STORAGE.

Traditionally users have used adhesive paper labels on their storage racks to provide sample traceability. However over a period of time, subject to repeated freeze thaw cycles and organic solvents which samples are often dissolved in, paper labels commonly degrade or fall off.

Each white Loborack-96 or Roborack-96 is laser etched with a guaranteed unique 1D barcode ensuring absolute sample traceability. The white rack with laser etched bar code provides an excellent bar coding readability and minimises the risk of incorrect sample identification.

The robust designed laser etched Loborack-96 or Roborack-96 incorporates a wealth of novel features enabling turn-key integration into almost any laboratory automation system.

Manufactured from hi-grade solvent resistant polypropylene,

both racks meet the globally recognized Society of Biomolecular Screening (SBS) standard footprint ensuring compatibility with most automated sample handling systems.


Available features:

- Laser etched unique barcode
- Sterile, gamma radiated
- Only white colour



INFORMATION

- Laser etched unique barcode
- Industry standard Society of Biomolecular Screening (SBS) footprint
- Secure and safe locking system
- Absolute consistency
- Working temperature ultra low
- Automation friendly
- Resists high chemical solvents
- Class 7 clean-room production
- RNase/DNase and Pyrogene free

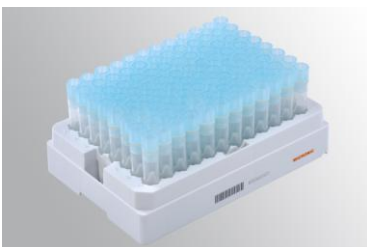


MICRONIC | SAMPLE STORAGE | RACKS

SPECIFICATIONS



MATERIAL	HIGHEST PURITY (MEDICAL GRADE) POLYPROPYLENE
SBS	8X12 MICROPLATE STANDARD (9MM PITCH)
RESISTANCE	RESISTS HIGH CHEMICAL SOLVENTS
AUTOCLAVABLE	YES

TECHNICAL DRAWING

I) Roborack-96

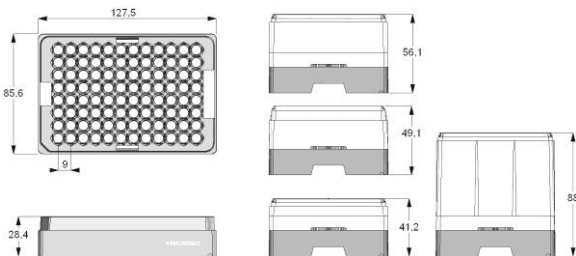

II) Lobarack-96

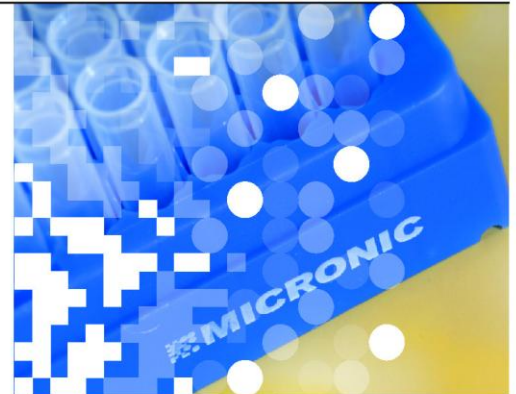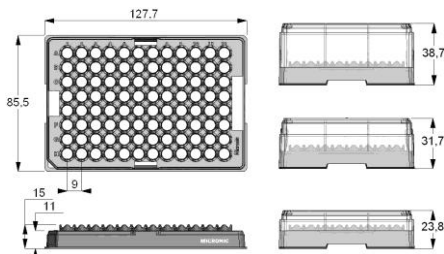

INFORMATION

- BC = Barcoded
- WS-side = West Short side

ORDERING INFORMATION

■ Laser Etched Roborack-96 or Lobarack-96		
MPW51001BC3	Roborack-96 White, standard cover (BC WS-side) for 1.40ml tubes with TPE or EVA pushcap	box/10
MPW51012BC3	Roborack-96 White with high cover (BC WS-side) for 1.40ml Screw Cap tubes	box/10
MPW51016BC3	Roborack-96 White with extra high cover (BC WS-side) for 2.50ml tubes with TPE or EVA pushcaps	box/10
MPW51004BC3	Lobarack-96 White barcoded with cover (BC WS-side) for 0.50ml tubes with screwcap or 0.75ml tubes TPE or EVA pushcaps	box/10
MPW51013BC3	Lobarack-96 White, with low cover (BC WS-side) for 0.50ml tubes with TPE Caps	box/10
MPW51015BC3	Lobarack-96 White, with high cover (BC WS-side) for 0.75ml tubes with Screwcap	box/10



Please visit our website www.Micronic.com or www.integra-biosciences.de

INTEGRA Biosciences GmbH
Ruhberg 4 • 35463 Fernwald
Tel. 06 404 / 809 / 0
Fax 06 404 / 809 / 251
info@integra-biosciences.de
www.integra-biosciences.de

