


MICRONIC | SAMPLE STORAGE | RACKS

LINE OF STORAGE RACKS-96

ARE YOU USING STORAGE TUBES IN 96-WELL CONFIGURATION? IF SO, THEN MICRONIC OFFERS A VARIETY OF STORAGE RACK SOLUTIONS. DEPENDING ON YOUR APPLICATION THE FOLLOWING STORAGE RACKS ARE AVAILABLE:

- ROBORACK-96
- LOBORACK-96
- COMORACK-96 OR COMORACK-24

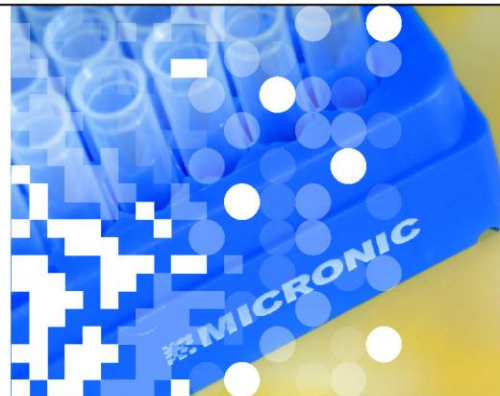


Roborack-96



Differences

- PP = Polypropylene
- Automation friendly
- Working temperature -80°C to +80°C
- SBS format 8x12 microplate standard (9mm pitch)
- For 1.10ml, 1.40ml and 2.50ml storage tubes
- With transparent PP lockable cover
- Special caps for Lyophilisation
- Autoclavable

Loborack-96

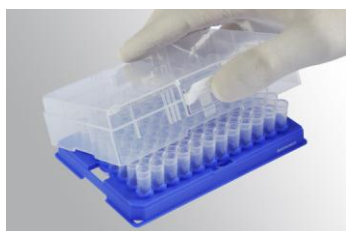

Differences

- PP = Polypropylene
- Automation friendly
- Working temperature -80°C to +80°C
- SBS format 8x12 plate standard (9mm pitch)
- For 0.50ml and 0.75ml storage tubes
- With transparent PP lockable cover
- Special caps for Lyophilisation
- Autoclavable

Comorack-96 or Comorack-24



Differences

- PC = Polycarbonate
- Robust designed and durable rack, can be re-used
- Working temperature -80°C to +80°C
- SBS format 8x12 microplate standard (9mm pitch) and 4x6 microplate standard
- For 1.10ml, 1.40ml and 7.50ml storage tubes
- With transparent PC lockable cover
- Autoclavable, 100 times at 121 °C for 15 minutes

INFORMATION

- Industry standard Society of Biomolecular Screening (SBS) footprint
- Class 7 clean room production
- RNase, DNase and Pyrogene free
- High chemical resistance
- Automation-friendly




MICRONIC | SAMPLE STORAGE | RACKS

EXPLANATION DIFFERENCES

Cover lock system, easy locking and unlocking mechanism for single-hand use



Sure Shot wells, automation-friendly, tapered wells allow precise tube entry into racks

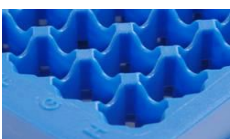

Twistlock, notches in well bottom prevent tubes from turning during screw-capping



Open rack-bottom, openings around tube wells decrease time of freeze and thaw cycles



Black racks, complete line of storage racks to keep your light sensitive samples safe.



Snap tubes, tube-locking mechanism holds tubes safely in the rack



Barcode Service, all Micronic storage racks are available bar-coded and sterile



Conductive racks: All Micronic polypropylene racks are available conductive for improved liquid-level sensing



Please visit our website www.Micronic.com or www.integra-biosciences.de

