

MICRONIC | SAMPLE STORAGE | RACKS

COMORACK-24

DO YOU WANT AN EFFECTIVE AND RELIABLE LOGISTICS AND STORAGE SYSTEM FOR YOUR VALUABLE SAMPLES? IF SO, THEN MICRONIC OFFERS YOU A ROBUST STORAGE RACK. THE COMORACK-24 HAS PROVEN TO BE A RELIABLE AND DURABLE RACK FOR MORE THAN 20 YEARS.

The Comorack-24 can hold 24 individual tubes (4x6), providing for the individual retrieval of your high valuable samples. This storage rack is specially designed for the Micronic 7.50ml sample storage tubes.

The Comorack-24 is based on the global recognised standard footprint from the Society of Bimolecular Screening (SBS).

The dimensions of the Comorack-24 are slightly different, as this rack was developed even before the SBS became the standard.

Manufactured to the industry leading strict tolerances, the precision and rack-to-rack consistency of Micronic products maximises operational uptime when used with automated rack handling systems. The rack is manufactured of the highest-purity polycarbonate, which resists many organic solvents.

This storage rack can be autoclaved up to 100 times. The Comorack-24 can be closed securely with a transparent lid. With the unique and easy-to-use orange locking system the lid stays perfectly closed throughout long-term storage and transportation.

Available features:

- Pre-labelled barcode
- Sterile, gamma radiated



INFORMATION

- Autoclavable, up to 100 times
- Secure and safe locking system
- Absolute consistency
- Ultra low Storage temperature
- Available with custom specific barcodes
- RNase & DNase and pyrogene free
- Class 7 clean room production




MICRONIC | SAMPLE STORAGE | RACKS

SPECIFICATIONS

MATERIAL	HIGHEST PURITY (MEDICAL GRADE) POLYCARBONATE
SBS	4X6 MICROPLATE STANDARD
RESISTANCE	RESISTS HIGH CHEMICAL SOLVENTS
AUTOCLAVABLE	YES



TECHNICAL DRAWING

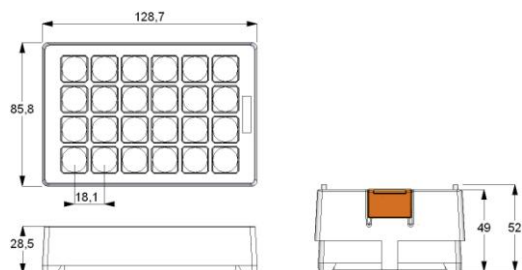

ORDERING INFORMATION

Comorack-24

MP23100	PC Comorack-24 with Cover	10 racks
MP23100BC	PC Comorack-24 with Cover (barcoded)	10 racks



Please visit our website www.Micronic.com or www.integra-biosciences.de



INTEGRA Biosciences GmbH
 Ruhberg 4 • 35463 Fernwald
 Tel. 06 404 / 809 / 0
 Fax 06 404 / 809 / 251
 info@integra-biosciences.de
 www.integra-biosciences.de

