

MICRONIC | SAMPLE STORAGE | STARTERPACKS

STARTERPACK WITH SCREW CAPS

DO YOU WANT TO GET STARTED USING 2D CODED SCREW CAP TUBES TO TRACK AND STORE YOUR VALUABLE SAMPLES? IF SO, THEN MICRONIC OFFERS YOU THE POSSIBILITY TO START USING A FULLY TRACEABLE STORAGE SOLUTION AT A LOW INVESTMENT. BUYING THE STARTERPACK CAN SAVE YOU UP TO 25% COMPARED TO BUYING THE PRODUCTS INDIVIDUALLY.

The Micronic Starterpack with Screw Caps is a complete package that ensures you a safe and traceable lab-process for your valuable samples. It contains everything needed to start using 2D coded sample storage Screw Cap tubes, enabling laboratory workers to ensure a secure sample logistics system and eliminate the costly possibility of false sample identities.

Micronic has Screw Cap Starterpacks available for tubes from volumes of 0.50ml, 0.75ml, 1.10ml and 1.40ml. The tubes are available with an inner V- or U-bottom. Each Starterpack includes 40 racks with 2D coded tubes (Tracker Data-Matrix or TraXis) in Lobarack-96 or Roborack-96.

The Single Rack Scanner (mvs) is the ideal solution for manual reading racks filled with 2D coded (Tracker Data-Matrix or TraXis) tubes within 15 seconds. The accompanying software is easy to use and gives a reliable decoding of the 2D codes on the tubes.

For a secure sealing of the tubes, the Starterpack contains enough Screw Caps for capping all 40 racks with tubes. The manual Screw Cap tool helps you with the capping of your tubes.

When individual tubes need to be picked from the rack, the handy Tube Selector is a perfect solution.


Available features

- Available sterile
- Barcode service
- Snap Tubes, locks tubes into the wells in the rack to prevent sample loss from overturned racks
- Pre-capped, with Screw Caps



INFORMATION

- Save up to 25%
- Enables users to start storing and sealing up to 4,000 samples
- Available in volumes of 0.5ml, 0.75ml, 1.10ml or 1.40ml
- Single Rack Scanner reads a full rack in 15 seconds
- 2D Decoding software included (Tracker Data-Matrix or TraXis)
- Class 7 clean room production
- RNase/DNase and Pyrogene free

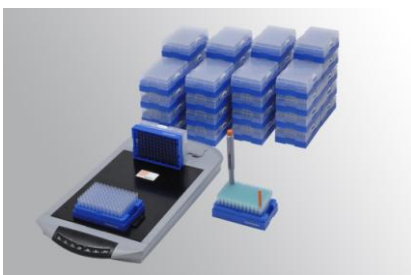

MICRONIC | SAMPLE STORAGE | STARTERPACKS

MINIMUM COMPUTER REQUIREMENTS

PROCESSOR	PENTIUM II 200MHZ
INTERNAL MEMORY	128 MB RAM, 2 GB HARD DISK
DRIVE	CD ROM
OPERATING SYSTEM	WINDOWS 2000, SP 2 OR WINDOWS XP
USB PORTS	2X



ORDERING INFORMATION

0.50ml Screw Cap tubes V-bottom Starterpacks


MP50102	0.50ml 2D Tracker Data-Matrix Coded Screw Cap tubes Starterpack	unit
MP52325	0.50ml 2D Tracker Screw Cap tubes in Lobarack-96	40 racks
MP53106	Screw Caps Transparent Blue	3,840 caps
MP54004	Manual Screw Cap remover	10 pcs
MP55114	Single Rack Scanner (mvs)	1 unit
MP66000	Tube Selector (Manual tube picker)	1 unit
MP50002	0.50ml 2D TraXis Coded Screw Cap tubes Starterpack	unit

1.10ml Screw Cap tubes V-bottom Starterpacks

MP50103	1.10ml 2D Tracker Data-Matrix Coded Screw Cap tubes Starterpack	unit
MP50003	1.10ml 2D TraXis Coded Screw Cap tubes Starterpack	unit

1.40ml Screw Cap tubes V-bottom Starterpacks

MP50104	1.40ml 2D Tracker Data-Matrix Coded Screw Cap tubes Starterpack	unit
MP50004	1.40ml 2D TraXis Coded Screw Cap tubes Starterpack	unit

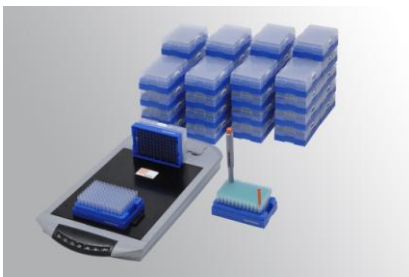

Please visit our website www.Micronic.com or www.integra-biosciences.de



INTEGRA Biosciences GmbH
 Ruhberg 4 • 35463 Fernwald
 Tel. 06 404 / 809 / 0
 Fax 06 404 / 809 / 251
info@integra-biosciences.de
www.integra-biosciences.de

