

MICRONIC | SAMPLE HANDLING | EQUIPMENT

MICRONIC MANUAL CAPMAT SEALER

DO YOU REQUIRE A SIMPLE AND COMPLETE SOLUTION TO CLOSE YOUR SAMPLES WITH TPE PLUG STYLE CAPS? THEN, THE MANUAL CAPMAT SEALER IS THE IDEAL MANUAL CAPPING SOLUTION FOR YOU. SAMPLE INTEGRITY IS A KEY REQUIREMENT FOR THE LONG-TERM STORAGE OF VALUABLE SAMPLES.

The Micronic Manual Capmat Sealer guarantees an equal sealing of plug style tubes in a 96-well configuration. This manual device requires no installation and is simple to operate.

The Manual Capmat Sealer can be used to seal any TPE cap on the Micronic Roborack-96, Lobarack-96 and Comorack-96. Using racks with an SBS-footprint minimises the storage space of your samples.

Standard for 1.40ml tubes in all configurations like: non-coded, Alphanumeric coded and 2D coded (Tracker Data-Matrix and TraXis). For different sizes in a 96-well configuration, adapter blocks are required for height adjustment. For the 0.50ml and 0.75ml tubes, we offer another type of Manual Capmat Sealer to guarantee easy capping.

The Micronic Capmats and Capclusters are recommended for the sealing of your valuable samples.

The Capmat exists out of a retaining film that holds 96 fixed TPE caps for sealing 96 tubes in a rack.



The Capcluster enables you to seal 96 tubes as easily as applying a Capmat, simply by removing the retaining film in order to leave each tube individually capped.

The functionality of the Capmat sealer is developed in conjunction with our users. Micronic has more than 25 years of experience in the life science market.



INFORMATION

- Standard for all 1.4ml tubes
- Easy installation
- Integrity of your sample ensured
- Robust design
- One slide capping method

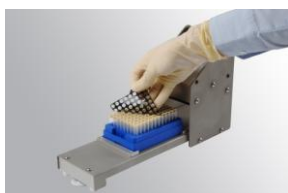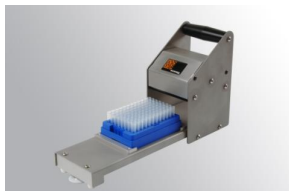


MICRONIC | SAMPLE HANDLING | EQUIPMENT

SPECIFICATIONS:

MATERIAL: STAINLESS STEEL
DIMENSIONS: 350CM L x 206 H x 114 W
WEIGHT: 6 KG

TECHNICAL DRAWING:

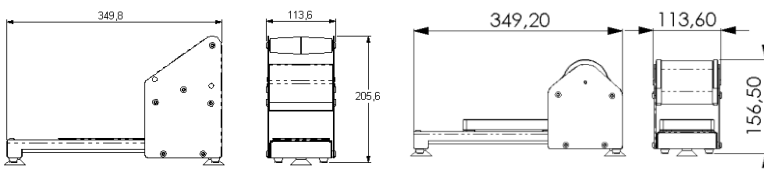

ORDER INFORMATION:


Manual Capmat Sealer

MP55005	Manual Capmat Sealer 1.40ml and 1.10ml Tubes for EVA & TPE	unit
MP55006	Manual Capmat Sealer 0.50 and 0.75ml Tubes for EVA & TPE	unit


Adapters for Manual Capmat Sealer (ordered separately)

MP55702	Manual Capmat Sealer Adapter for 0.50ml Tubes	unit
MP55700	Manual Capmat Sealer Adapter for 0.75ml Tubes	unit
MP55701	Manual Capmat Sealer Adapter for 1.10ml Tubes	unit

MP55006

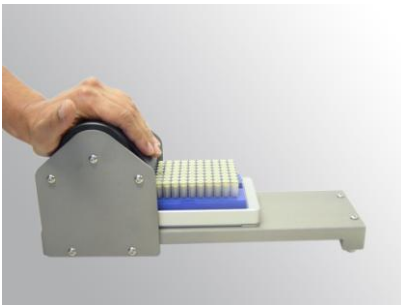

MP55005

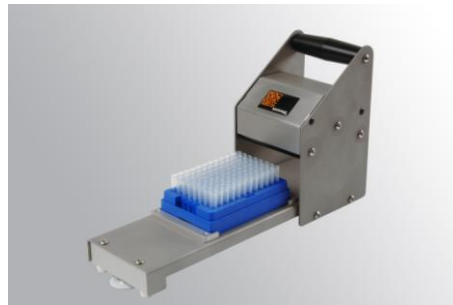

Please visit our website www.Micronic.com or www.integra-biosciences.de



INTEGRA Biosciences GmbH
 Ruhberg 4 • 35463 Fernwald
 Tel. 06 404 / 809 / 0
 Fax 06 404 / 809 / 251
 info@integra-biosciences.de
 www.integra-biosciences.de

