

MICRONIC | SAMPLE HANDLING | EQUIPMENT

ELECTRIC CAPPER

MAKE YOUR CAPPING EASY. MICRONIC OFFERS AN ELECTRIC CAPPER FOR THE SEALING OF A VARIETY OF TUBE SIZES WITH CAPMATS, CAPBANDS AND CAPCLUSTERS. THE ELECTRIC CAPPER IS A LOW COST DEVICE WITH ALL THE ROBUSTNESS AND EASY-TO-USE FEATURES THAT YOU CAN EXPECT FROM A BENCHTOP DEVICE.

Built to meet the special individual requirements of our customers, this Electric Capper can cap all Micronic plug style caps: capmats, capbands and capclusters made of EVA or TPE on Micronic tubes and racks. The capper is also compatible with other brands of tubes and caps (please contact us for a list of compatible brands).

The Electric Capper has the robustness and easy-to-use features that you expect from a benchtop device. For capping, simply place the rack filled with tubes in the drawer, and then place the capmat, capband or capcluster on the tubes, close the drawer and press the button.

The controlled pressure system will press all the caps down evenly to ensure the integer capping of your tubes. Simply pull the drawer out and retrieve the rack.

The Electric Capper has the flexibility to cap both full and partially-full racks of tubes. This is possible due to the controlled pressure system that caps each tube individually.



INFORMATION

- Caps the entire rack of 96 tubes in only 10 seconds
- Compatible with the complete Micronic range of tubes and racks
- Also compatible with other brands of tubes and caps (please contact us for a list of compatible brands)
- Works standard for 1.40ml tubes, but smaller sized adapters are available
- Works on full and partially-full racks



MICRONIC | SAMPLE HANDLING | EQUIPMENT

SPECIFICATIONS

COMPATIBILITY	ALL SBS STANDARD RACKS
SPEED	10 SECONDS PER RACK
DIMENSIONS	38CM X 15CM X 33CM
WEIGHT	13KG
VOLTAGE	115 - 230 V




ORDERING INFORMATION

Electric Capper

MP35015	Electric Capper	unit
---------	-----------------	------

Adapters for Electric Capper

MP55702	Adapter for 0.50ml Micronic tubes	box/1
MP55700	Adapter for 0.75ml Micronic tubes	box/1
MP55701	Adapter for 1.10ml Micronic tubes	box/1

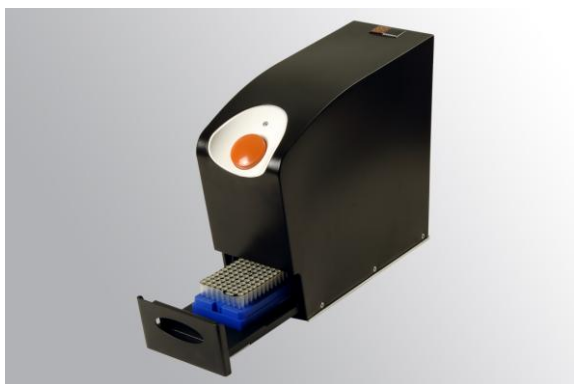

Please visit our website www.Micronic.com or www.integra-biosciences.de



INTEGRA Biosciences GmbH
 Ruhberg 4 • 35463 Fernwald
 Tel. 06 404 / 809 / 0
 Fax 06 404 / 809 / 251
info@integra-biosciences.de
www.integra-biosciences.de

