

MICRONIC | SAMPLE HANDLING | EQUIPMENT

TRACK-IT DATABASE SOFTWARE

DO YOU REQUIRE KEEPING TRACK OF YOUR SAMPLE STORAGE TUBES? TRACK-IT CAN HELP YOU MANAGE YOUR SAMPLES STORAGE AND KEEP TRACK OF ALL ACTIVITIES AND ALIQUOTS. IT WILL HELP YOUR LAB TO CARRY OUT BIOBANKING QUICKLY AND EASILY. NO MORE PAPER RECORDS OR LOSE EXCEL FILES.

Track-IT is a new Laboratory Sample Management System (LSMS) that allows you to organise and track individual samples in your library. Track-IT is compatible with barcode scanning, in turn allowing you to catalogue samples that have been scanned by lab equipment, such as readers and scanners.

The database is customisable to match the workflow of your laboratory, as well as to match the unique aspects of each of your samples.


The main concept behind Track-IT is that each scientist in the lab can manage data and information quickly and make it available to the rest of the lab community.

Via Intranet support, data can be accessed from any connected computer in the lab, which accelerates information retrieval.

Laboratory staff in turn spends less time on recurrent and time-consuming operations and dedicate more precious time to research activity.

Features and benefits:

- Central information database
- Quick information retrieval
- Alerts system
- Cross-database searches
- Uniform data presentation
- Perpetuation of information & Traceability
- Data reports export (screen and MS-Excel)
- Automatic data imports
- Versatile administration
- Upgradable
- 1 licence per database covers all the lab users

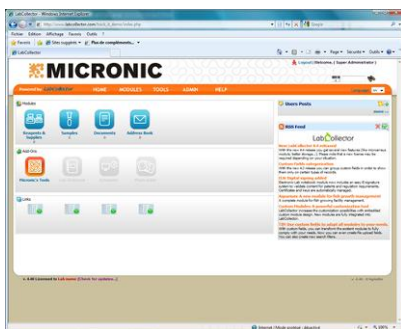

MICRONIC | SAMPLE HANDLING | EQUIPMENT

SPECIFICATIONS



ALL-IN-ONE AUTOMATIC WINDOWS INSTALLER


SERVER SYSTEM REQUIREMENTS: WINDOWS 2000/XP/VISTA, 2X USB PORTS


CLIENT SYSTEM REQUIREMENTS: WINDOWS/MACOS X/LINUX + INTERNET BROWSER (INTERNET EXPLORER 6+, FIREFOX 2+, SAFARI 2+, OPERA 8+)

Track-IT installs on one computer that becomes the server. All other users can connect to it from any Internet browser.



ORDERING INFORMATION

	Track-IT Database Software	
RH81001	Track-IT Laboratory Sample Management System (LSMS)	unit

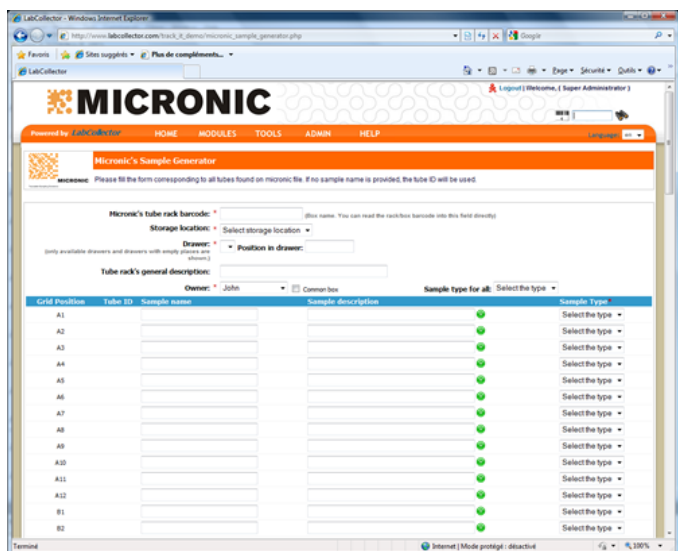

Please visit our website www.Micronic.com or www.integra-biosciences.de



INTEGRA Biosciences GmbH
 Ruhberg 4 • 35463 Fernwald
 Tel. 06 404 / 809 / 0
 Fax 06 404 / 809 / 251
info@integra-biosciences.de
www.integra-biosciences.de

